




bright System

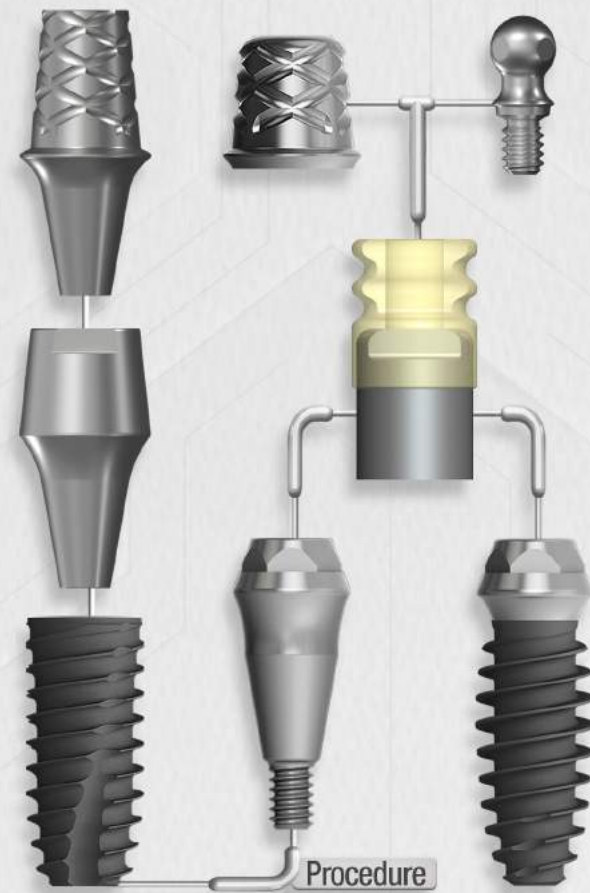

간단하다

단순한 술식, 간단한 보철

The **Bright**

Manufacturer **FIMS Co., Ltd.**
56, Changnyong-daero 256beon-gil, Yeongtong-gu, Suwon-si, Gyeonggi-do

Distributed by **The Bright Co., Ltd**
F10, 21, Teheran-ro 87-gil, Gangnam-gu, Seoul

bright

Bone Level Tissue Level



www.brightz.co.kr

bright System

bright System



www.brightz.co.kr

bright Implant

Index

Fixture Line-up

Fixture (Bone, Tissue Level)	12
-------------------------------------	----

Prosthetic Line-up

Scan Abutment	16
Metal Link Abutment	17
Temporary Abutment	18
Screw Abutment	19
Angled Screw Abutment	20
Screw Abutment Components	21
Screw Scan Comfort Cap	21
Screw Pick-up Impression Coping Screw	21
Screw Metal Link	22
Screw PEEK Link	22
Temporary Cylinder TI	22
Temporary Cylinder PEEK	23

Overdenture Line-up


Screw Ball Cap	24
Screw Ball Cap Socket Spacer & O-ring	24
Screw Ball Impression Coping	25
Screw Ball Analog	25
Screw Ball Cap Socket Spacer	25
Lab Analog	25


Kit & Instrument Line-up

Drilling Sequence	27
--------------------------	----


Kit Components	
Kit Includes Instrument	



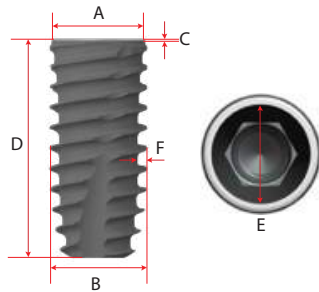





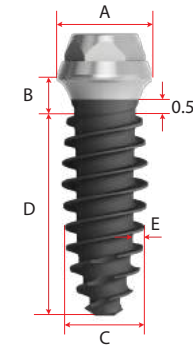

Bone Level Fixture



Fixture (Mount Free)							
A	Platform Diameter(mm)	Ø3.65	Ø3.65	Ø3.80	Ø4.00	Ø4.50	Ø5.90
B	Body Diameter(mm)	Ø2.90	Ø3.50	Ø4.00	Ø4.40	Ø5.00	Ø5.00
C	Bevel Height (mm)	2.00					
	L: 7 / 9 / 11 / 13	0	0	0.05	0.10	0.25	0.65
D	Total Length(mm)	7 / 9 / 11 / 13					
E	Abutment Interface(mm)	Ø3.30	Ø3.30	Ø3.30	Ø3.30	Ø3.30	Ø3.30
F	Thread Depth(mm)	0.30	0.35	0.40	0.45	0.50	0.50
Cap Color							
		Gray	Yellow	Green	Blue	Red	Orange
Selection Guideline		Anterior	Anterior	Premolar	Premolar	Molar	Molar

Tissue Level Fixture



Fixture (Mount Free)						
A	Platform Diameter(mm)	Ø3.40	Ø3.40	Ø4.20	Ø4.20	Ø4.20
B	Gingival Height(mm)	1.50	1.50	1.50	1.50	1.50
C	Body Diameter(mm)	Ø2.50	Ø3.00	Ø3.50	Ø4.00	Ø5.10
D	Total Length(mm)	7 / 9 / 11 / 13				
E	Thread Depth(mm)	0.30	0.45	0.50	0.55	0.65
Cap Color						
		Gray	Gray	Yellow	Green	Red
Selection Guideline		Anterior	Anterior	Premolar	Molar	Molar

bright Implant

Future of implants



SIMPLE

Simplify Surgical
Procedure



EASY

Easier Handling With
the Abutment



STRONG

Strength &
Joint Stability

S.L.A. Surface

PREDICTABLE

High Predictability &
Proven Surface Treatment

Simplified drill protocol



bright Bone / Tissue level
Common KIT



Digital Impression with
Healing Abutment

Repair 적은 임플란트

미래를 보는 임플란트

쉽고 간단하게 수술하세요.

걱정없이 안전하게 수술하세요.

bright Line-up

제품구매 및 bright 임상 증례
정보를 확인할 수 있습니다.




Line-up




Bright (Bone, Tissue Level) System

Prosthetic Line-up

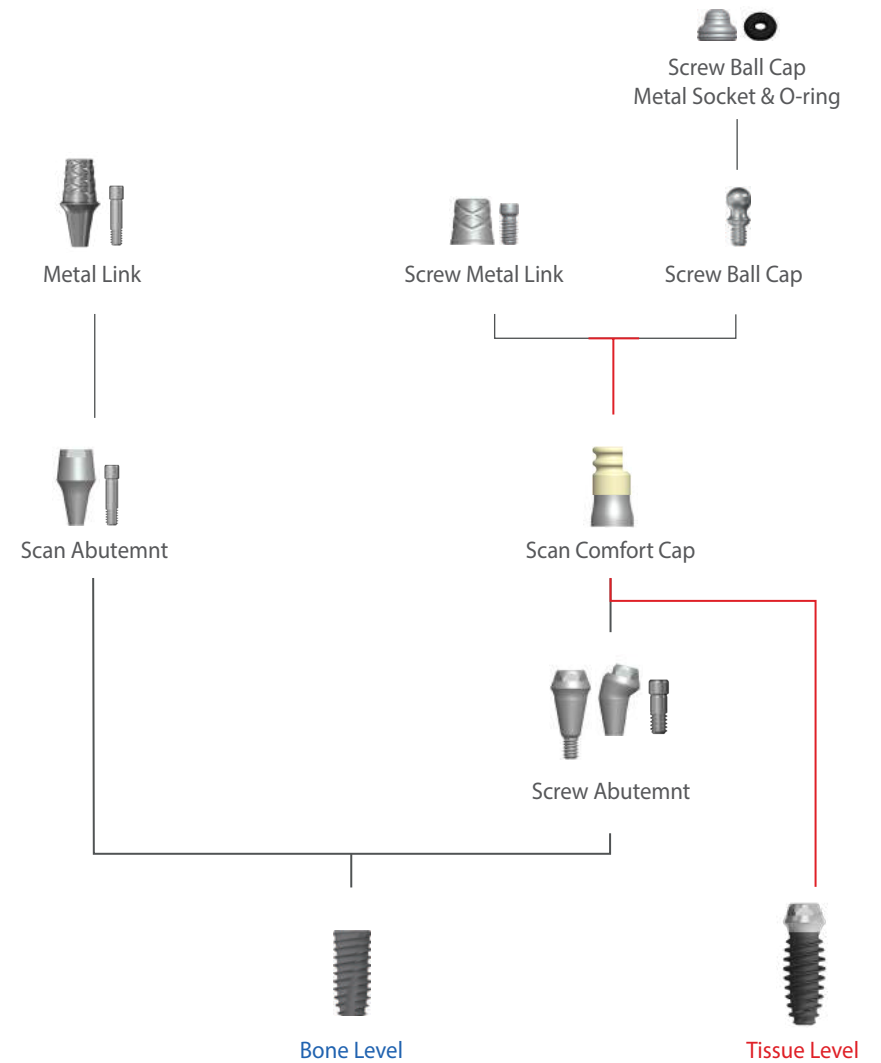


Bone Level Fixture

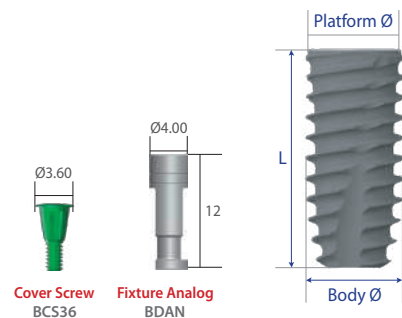


※ Hex driver: Use no more than 5N·cm of torque when screwing a Cover Screw to a fixture.

• 10°, 3.70mm Deep Connection

No rotation even in non-hex after adaptation of abutment.

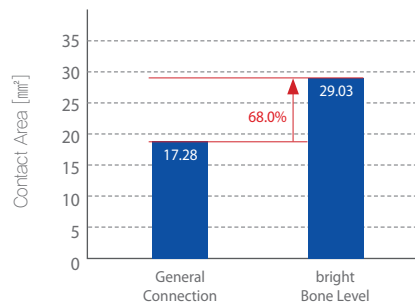

• Less Sink-Down

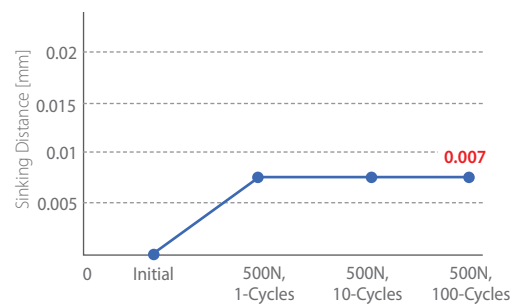

*당사 시험 기준

BFX30

Platform Ø3.65

Body Ø2.90

L	7.0	7.0	9.0	11.0	13.0
Art.No	BFX 30 07 B	BFX 30 07	BFX 30 09	BFX 30 11	BFX 30 13

BFX35

Platform Ø3.65

Body Ø3.50

L	7.0	7.0	9.0	11.0	13.0
Art.No	BFX 35 07 B	BFX 35 07	BFX 35 09	BFX 35 11	BFX 35 13

BFX40

Platform Ø3.80

Body Ø4.00

L	7.0	7.0	9.0	11.0	13.0
Art.No	BFX 40 07 B	BFX 40 07	BFX 40 09	BFX 40 11	BFX 40 13

BFX45

Platform Ø4.00

Body Ø4.40

L	7.0	7.0	9.0	11.0	13.0
Art.No	BFX 45 07 B	BFX 45 07	BFX 45 09	BFX 45 11	BFX 45 13

BFX50

Platform Ø4.50

Body Ø5.00

L	7.0	7.0	9.0	11.0	13.0
Art.No	BFX 50 07 B	BFX 50 07	BFX 50 09	BFX 50 11	BFX 50 13

BFX60

Platform Ø5.90



Body Ø5.05

L	7.0	7.0	9.0	11.0	13.0
Art.No	BFX 60 07 B	BFX 60 07	BFX 60 09	BFX 60 11	BFX 60 13



Tissue Level

Fixture (Included Cover Screw)

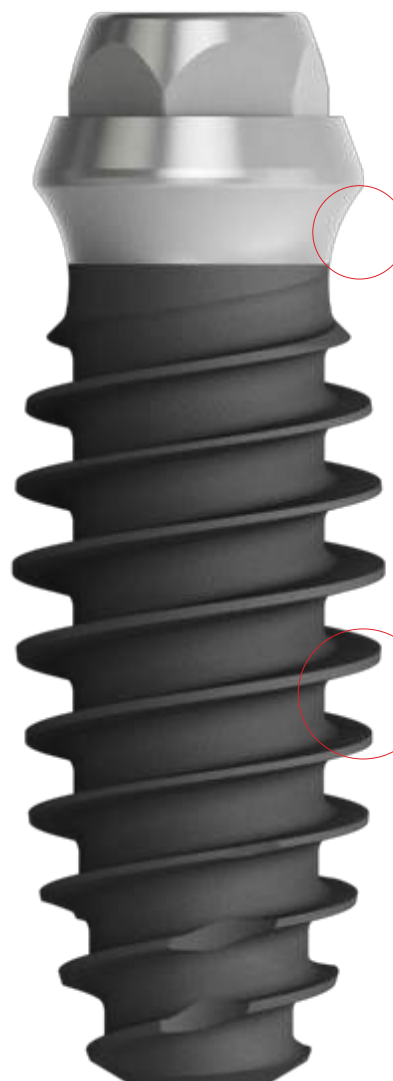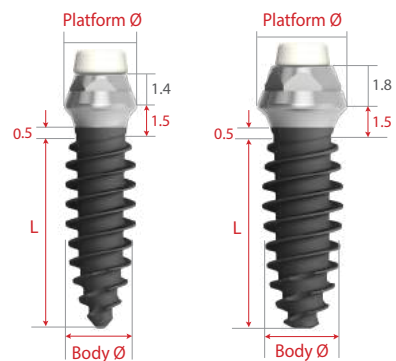

Hygiene zone

Minimal post height and adjustable placement depth according to soft tissue thickness



• 2-Platforms

For All Indications, prosthetic profile

- Platform Ø3.4 for Anterior
- Platform Ø4.2 for Posterior

Thread

Optimal Primary Stability and Bone-Implant Contact Area

NTI25

Platform Ø3.4

Body Ø2.5

*G/H 1.5

L	7.0	7.0	9.0	11.0	13.0
Art.No	NTI 34 25 07 B	NTI 34 25 07	NTI 34 25 09	NTI 34 25 11	NTI 34 25 13

NTI30

Platform Ø3.4

Body Ø3.1

*G/H 1.5

L	7.0	7.0	9.0	11.0	13.0
Art.No	NTI 34 30 07 B	NTI 34 30 07	NTI 34 30 09	NTI 34 30 11	NTI 34 30 13

NTI35

Platform Ø4.2

Body Ø3.5

*G/H 1.5

L	7.0	7.0	9.0	11.0	13.0
Art.No	NTI 42 35 07 B	NTI 42 35 07	NTI 42 35 09	NTI 42 35 11	NTI 42 35 13

NTI40

Platform Ø4.2

Body Ø4.0

*G/H 1.5

L	7.0	7.0	9.0	11.0	13.0
Art.No	NTI 42 40 07 B	NTI 42 40 07	NTI 42 40 09	NTI 42 40 11	NTI 42 40 13

NTI50

Platform Ø4.2

Body Ø5.1


*G/H 1.5

L	7.0	7.0	9.0	11.0	13.0
Art.No	NTI 42 50 07 B	NTI 42 50 07	NTI 42 50 09	NTI 42 50 11	NTI 42 50 13

※ Note: The Cover screw must be loosened when connecting prosthesis
Loosen the cover screw counterclockwise using Hex driver.

Scan Abutment (Non-Hex)

Healing, Scan Jig

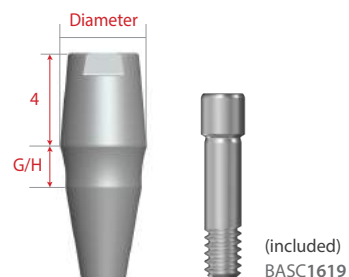

Diameter	Ø3.8		
G/H	1.0	2.0	3.0
Art.No	BIH38 10 50TN	BIH38 20 60TN	BIH38 30 70TN



Diameter	Ø5.0		
G/H	1.0	2.0	3.0
Art.No	BIH50 10 50TN	BIH50 20 60TN	BIH50 30 70TN




Diameter	Ø6.2		
G/H	1.0	2.0	3.0
Art.No	BIH62 10 50TN	BIH62 20 60TN	BIH62 30 70TN



※ Note: It is recommended to keep the torque level at 10 N-cm to tighten Abutment with fixture.

Metal Link Abutment (Non-Hex)

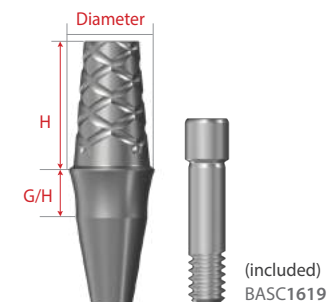


Diameter	Ø3.8		Ø5.0	
G/H	1.5	2.5	1.5	2.5
H	6.0	6.0	5.5	5.5
Art.No	BM38 15 N	BM38 25 N	BM50 15 N	BM50 25 N




※ Note: It is recommended to keep the torque level at 20~25 N-cm to tighten Abutment with fixture.

Abutment Screw


Art.No BASC1619

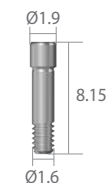

Temporary Abutment (Non-Hex)

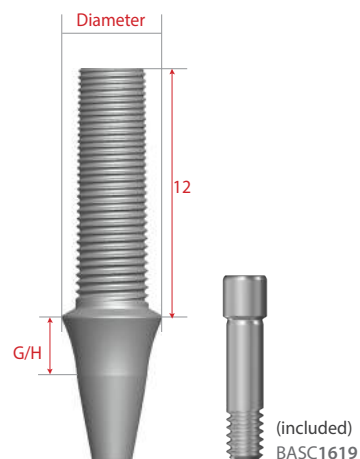


	Temporary Abutment TI		Temporary Abutment PEEK	
Diameter	Ø4.5			
G/H	1.0	3.0	1.0	3.0
Art.No	BTA45 10TN	BTA45 30TN	BTA45 10PN	BTA45 30PN

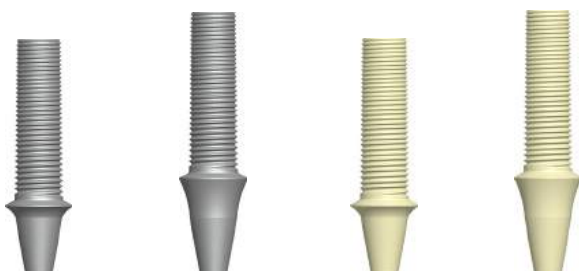

※ Note: It is recommended to keep the torque level at 15 N-cm to tighten Abutment with fixture.

Screw Abutment (Non-Hex)

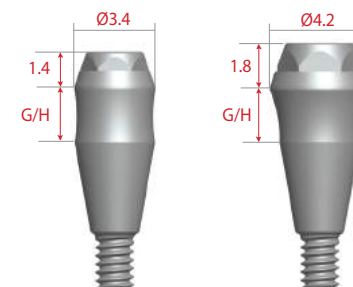


Diameter	$\varnothing 3.4$			
G/H	0.5	1.5	2.5	3.5
Art.No	BSA 34 05	BSA 34 15	BSA 34 25	BSA 34 35

Diameter	$\varnothing 4.2$			
G/H	0.5	1.5	2.5	3.5
Art.No	BSA 42 05	BSA 42 15	BSA 42 25	BSA 42 35

※ Note: It is recommended to keep the torque level at 20~25 N-cm to tighten Abutment with fixture.

Angled Screw Abutment (Non-Hex)

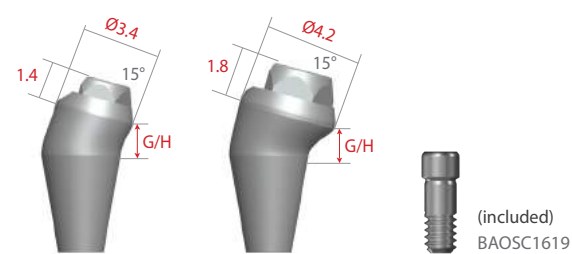

· Angled 15°


Diameter	Ø3.4		
G/H	0.5	1.5	2.5
Art.No	BASA34 15 05N	BASA34 15 15N	BASA34 15 25N

Diameter	Ø4.2		
G/H	0.5	1.5	2.5
Art.No	BASA42 15 05N	BASA42 15 15N	BASA42 15 25N

※ Note: It is recommended to keep the torque level at 20~25 N-cm to tighten Abutment with fixture.

Angled Screw Abutment Screw


Art.No BAOSC1619

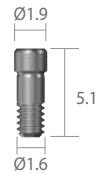

Screw Abutment Components (Hex)

Scan Jig, Impression Coping, Healing

* Include Comfort Cap Holder and Screw

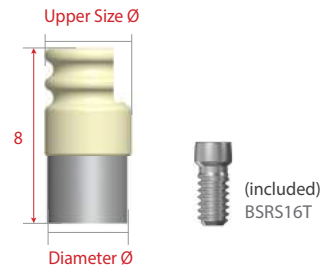

Screw Scan Comfort Cap


Application	Ø3.4		Ø4.2	
Diameter	Ø3.6		Ø4.5	
Upper Size Ø	3.2	4.0	4.0	5.0
Art.No	BISC34 27 HT	BISC34 36 HT	BISC42 36 HT	BISC42 48 HT

※ Note: It is recommended to keep the torque level at 10 N-cm to tighten Abutment with fixture.

Screw Pick-up Impression Coping Screw

Diameter	Ø2.5
Art.No	BSPS15

Screw Abutment Components (Hex)

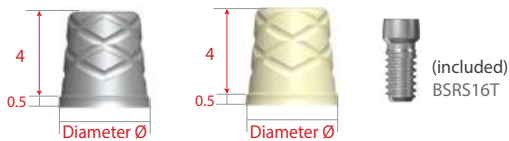

Screw Metal Link


Application	Ø3.4		Ø4.2	
Diameter	Ø3.4	Ø4.2	Ø4.2	Ø5.0
Art.No	BSM 34 34 05HT	BSM 34 42 05HT	BSM 42 42 05HT	BSM 42 50 05HT


Screw PEEK Link

Application	Ø3.4		Ø4.2	
Diameter	Ø3.4	Ø4.2	Ø4.2	Ø5.0
Art.No	BSM 34 34 05HP	BSM 34 42 05HP	BSM 42 42 05HP	BSM 42 50 05HP

※ Note: It is recommended to keep the torque level at 20~25 N-cm to tighten Abutment with fixture.

Screw Abutment Components (Hex)

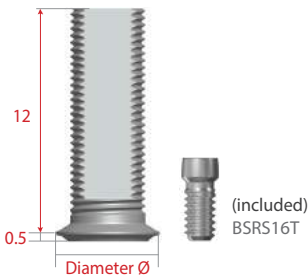


Temporary Cylinder TI		Temporary Cylinder PEEK		
Application	Ø3.4	Ø4.2	Ø3.4	Ø4.2
Diameter	Ø4.2	Ø5.0	Ø4.2	Ø5.0
Art.No	BSTC34 42 05HT	BSTC42 50 05HT	BSTC34 42 05HP	BSTC42 50 05HP
				

※ Note: It is recommended to keep the torque level at 15 N-cm to tighten Abutment with fixture.

Ti-Retaining Screw

Art.No BSR516T

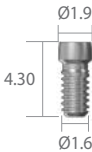

Screw Abutment Components


Overdenture

Screw Ball Cap

Diameter	Ø2.5
Art.No	BSBC25

· Delivery Holder is included


※ Note: It is recommended to keep the torque level at 30 N-cm to tighten Abutment with fixture.

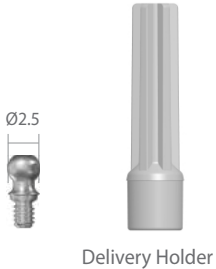

Screw Ball Cap Metal Socket & O-ring

Diameter	Ø5.5
Art.No	BSBMS

· Type 500~700gf

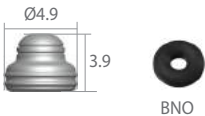

Screw Abutment Components

Overdenture

Screw Ball Impression Coping

Diameter	Ø4.05
L	4.2
Art.No	BICA



Screw Ball Analog

Diameter	Ø4.0
Art.No	BAL




Screw Ball Cap Socket Spacer

Diameter	Ø6.4
Art.No	BSKS




Lab Analog

Diameter	Ø3.4	Ø4.2
Art.No	BSDAN34	BSDAN42

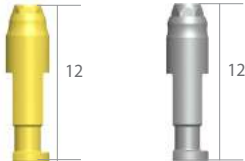

Kit & Instrument Line-up

Kit Components

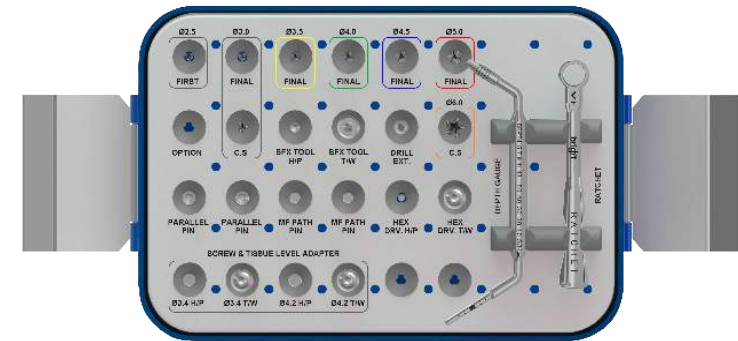

Bright Surgical Kit

Kit Includes

BSK

Guide Drill



BGD2433



BFD3033



BFD3533



BFD4033



BFD4533



BFD5033

Countersink



BCS3029

BFX TOOL H/P



BFFA28H

BFX TOOL T/W



BFD4033

Drill Extension



BDE

Countersink



BCS6027

Depth Gauge



BDG

Ratchet



BRCT

Parallel Pin



BPP164325

Path Pin



BEPP

Hex Driver



BHD25H



BHD25W

Screw & Tissue Level Adapter



BSA3421H



BSA3421W




BSA4221H



BSA4221W


Bone Level Drilling Sequence

Drilling Protocol

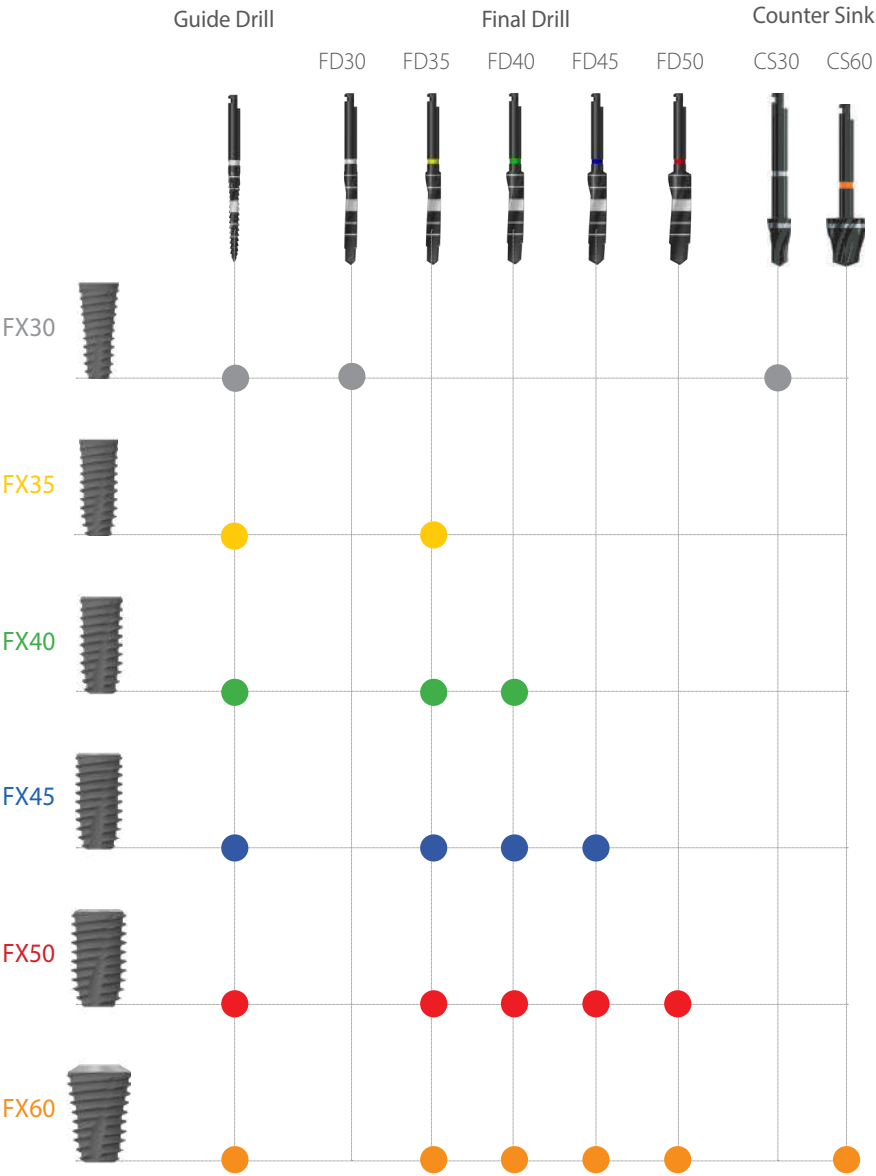

Tissue Level Drilling Sequence

Drilling Protocol

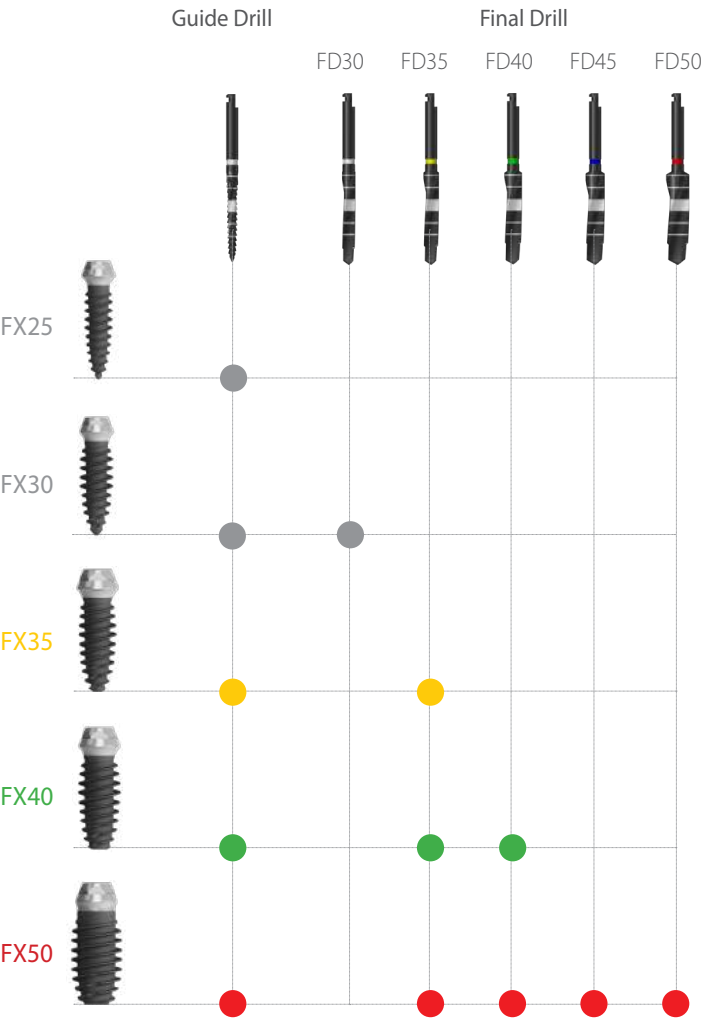

Instrument

First Guide Drill

Diameter	L	Art. No.
Ø 2.4	29.0	BGD 2529
	31.0	BGD 2531
	33.0	BGD 2533

Final Drill


Diameter	L	Art. No.
Ø 2.50	29.0	BFD 3029
	31.0	BFD 3031
	33.0	BFD 3033

Diameter	L	Art. No.
Ø 3.00	29.0	BFD 3529
	31.0	BFD 3531
	33.0	BFD 3533

Diameter	L	Art. No.
Ø 3.40	29.0	BFD 4029
	31.0	BFD 4031
	33.0	BFD 4033

Diameter	L	Art. No.
Ø 3.70	29.0	BFD 4529
	31.0	BFD 4531
	33.0	BFD 4533

Diameter	L	Art. No.
Ø 4.35	29.0	BFD 5029
	31.0	BFD 5031
	33.0	BFD 5033

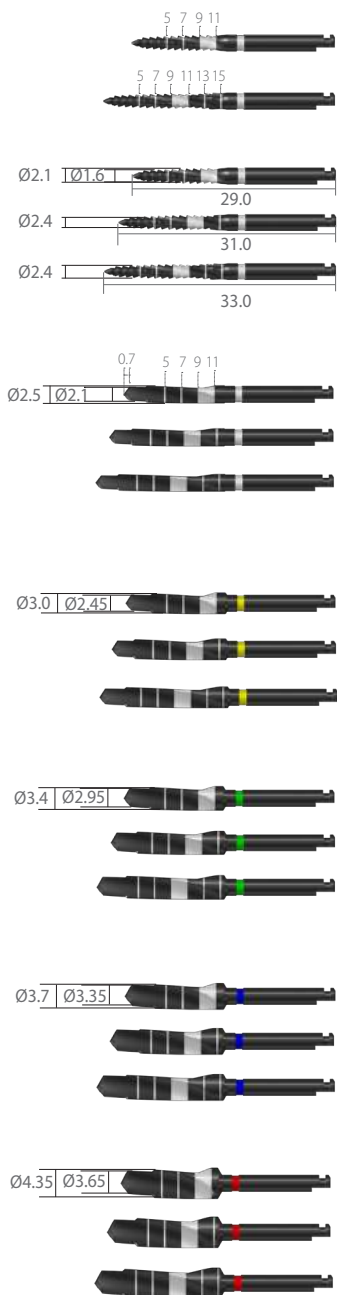


※ Note: Drill speed 1,000rpm, 30~45 N-cm with irrigation

Instrument

Countersink Drill


Diameter	L	Art. No.
Ø 3.65	29.0	BCS 3029
Ø 6.40	27.0	BCS 6027

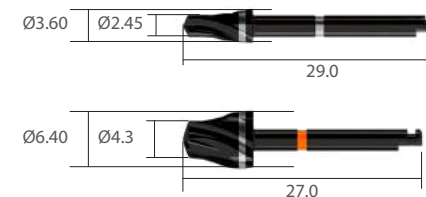

※ Note: Countersink Drill speed 20~60rpm, 30~45 N-cm with irrigation

Fixture Adapter

Type	L	Art. No.
Hand-Piece	27	BFA27H
	30	BFA30H
	32	BFA32H
Ratchet	24	BFA24W
	26	BFA26W
	29	BFA29W

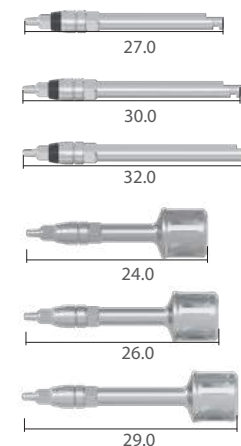

Parallel Pin

Diameter	Ø4.3
L	25
Art.No.	BPP164325




Path Pin

L	18.2
Art.No.	BFPP




Instrument


Hex Driver

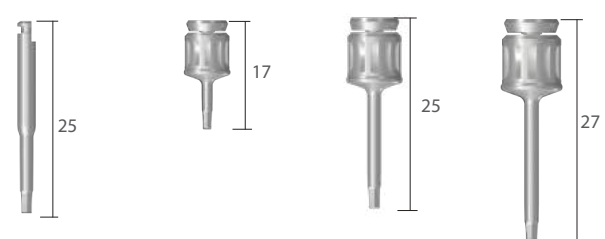Type	Hand-Piece	Ratchet		
L	25	17	25	27
Art.No.	BHD25H	BHD17W	BHD25W	BHD27W
				

Angled Hex Driver

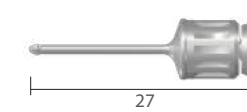Type	Ratchet
L	27
Art.No.	BAHD27W
	

Adapter for Screw Abutment

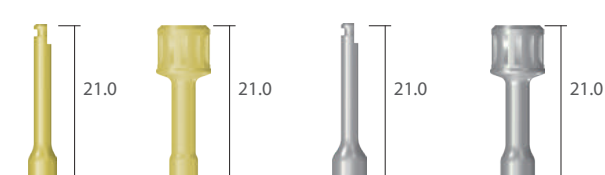Type	Hand-Piece	Ratchet	Hand-Piece	Ratchet
L	21	21	21	21
Art.No.	BSA3421H	BSA3421W	BSA4221H	BSA4221W
				

Instrument

Adapter for Ball Abutment

Type	Hand-Piece	Ratchet
L	21	21
Art.No.	BBA3821H	BBA3821W
		


Ratchet

Art.No. BRCT



Torque Wrench

Art.No. BTW

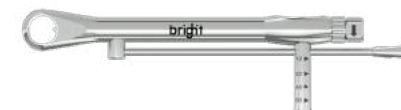

Depth Gauge

Art.No. BDG



bright System

bright System